

KitPipeP

KOO MEDICAL



Now with the new
ECO manometer



A **simple**
practice for a
great
benefit

KitPipeP

The complete range of kits and accessories, offers solutions both to hospital and home care therapeutical requirements. 8 kits including different mask sizes fit a range of patients from newborns to large adults. Standard cleaning, disinfection and sterilization methods **prevent infection transmission.**



Made of polycarbonate, the PipeP Kit features 9 resistors with different color-coded intensity for easier identification. **Two unidirectional silicon valves inside the valve body control inspiratory and expiratory flows** avoiding pressure alteration during the therapy. An arrow on the valve body **indicates the flow direction** during the breathing acts.

The valve mouthpiece can be either used **directly into the mouth**, or it **may be connected to masks of different sizes** (with a 22 mm inner diameter or 15 mm outer diameter connector). **Tracheostomy cannulas** with 15 mm outer diameter can be directly connected.



All Home Care Kits (KM-803HC, HCx) now with ECO manometer included for **PEP control from 0 to 20 cmH₂O.**

Indicazioni:

- **Secretion removal,**
- **Breathing volumes balance,**
- **Pulmonary rehabilitation following diseases, traumas or surgeries**

HomeCare Kit



Kit includes:

- Valve
- 9 resistors
- Mask
- T-connector with manometer connecting tube and ECO manometer
- Carrying bag

Hospital Kit



Kit includes:

- Valve
- 9 resistors
- T-connector with manometer connecting tube (manometer not included)

Resistors range




KM-803D	KM-803W	KM-803Y	KM-803B	KM-803G	KM-803O	KM-803M	KM-803C	KM-803P
1.5 mm	2.0 mm	2.5 mm	3.0 mm	3.5 mm	4.0 mm	4.5 mm	5.0 mm	6.0 mm
Black	White	Yellow	Blue	Green	Orange	Pink	Red	Purple


Codice	Descrizione	Pezzi per box
KM-803	Kit composition: valve, 9 resistors and T-connector with manometer connecting tube (manometer not included)	1 kit
KM-803HC	Kit composition: valve, 9 resistors, T-connector with manometer connecting tube, ECO manometer, mouthpiece and carrying bag	1 kit
KM-803HCI	Kit composition: valve, 9 resistors, T-connector with manometer connecting tube, ECO manometer, newborn 0 silicone mask and carrying bag	1 kit
KM-803HCO	Kit composition: valve, 9 resistors, T-connector with manometer connecting tube, ECO manometer, newborn 1 silicone mask and carrying bag	1 kit
KM-803HCN	Kit composition: valve, 9 resistors, T-connector with manometer connecting tube, ECO manometer, infant silicone mask and carrying bag	1 kit
KM-803HCP	Kit composition: valve, 9 resistors, T-connector with manometer connecting tube, ECO manometer, pediatric silicone mask and carrying bag	1 kit
KM-803HCA	Kit composition: valve, 9 resistors, T-connector with manometer connecting tube, ECO manometer, adult silicone mask and carrying bag	1 kit
KM-803HCL	Kit composition: valve, 9 resistors, T-connector with manometer connecting tube, ECO manometer, adult extra large silicone mask and carrying bag	1 kit
KM-1034	PVC Panda Mask, infant 0 – 3 years	50 PCS
KM-1033	PVC Panda Mask, pediatric 3 - 6 years	50 PCS
KM-8045	Silicone mask newborn 0	100 PCS
KM-8044	Silicone mask newborn 1	100 PCS
KM-8043	Silicone mask infant	100 PCS
KM-8042	Silicone mask pediatric	100 PCS
KM-8041	Silicone mask adult	100 PCS
KM-8040	Silicone mask extra-large	100 PCS
KM-803D	Resistors set 1,5 mm diameter, black	10 PCS
KM-803W	Resistors set 2,0 mm diameter, white	10 PCS
KM-803Y	Resistors set 2,5 mm diameter, yellow	10 PCS
KM-803B	Resistors set 3,0 mm diameter, blue	10 PCS
KM-803G	Resistors set 3,5 mm diameter, green	10 PCS
KM-803O	Resistors set 4,0 mm diameter, orange	10 PCS
KM-803M	Resistors set 4,5 mm diameter, pink	10 PCS
KM-803C	Resistors set 5,0 mm diameter, red	10 PCS
KM-803P	Resistors set 6,0 mm diameter, purple	10 PCS
KM-803K	Mixed resistors set, from 1,5 mm to 6,0 mm diameter	9 PCS
KM-8032	T-connector with 50 cm connecting tube and ECO manometer	50 PCS
KM-8035	T-connector with 50 cm connecting tube	50 PCS
KM-8051	Mouthpiece	50 PCS
57301537	Scaled manometer – 50 +50 cmH ₂ O	1 PC
57303499	Scaled manometer – 30 +30 cmH ₂ O	1 PC

Assembly instructions


Kit assembly

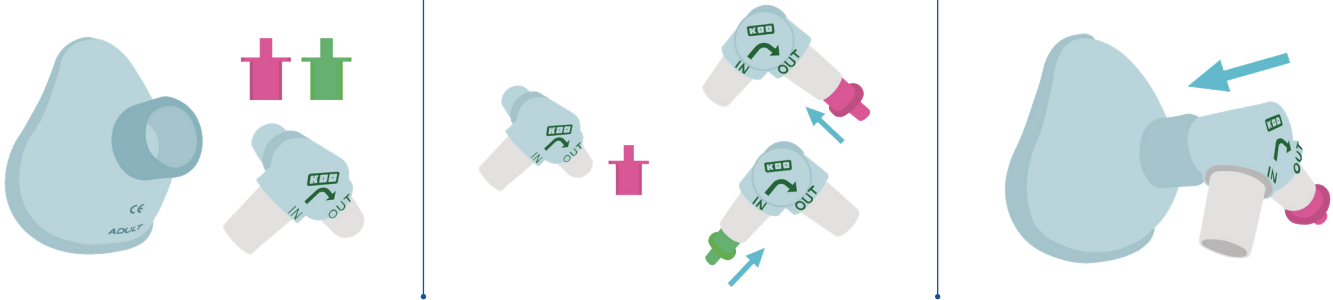

Kit assembly with standard manometer

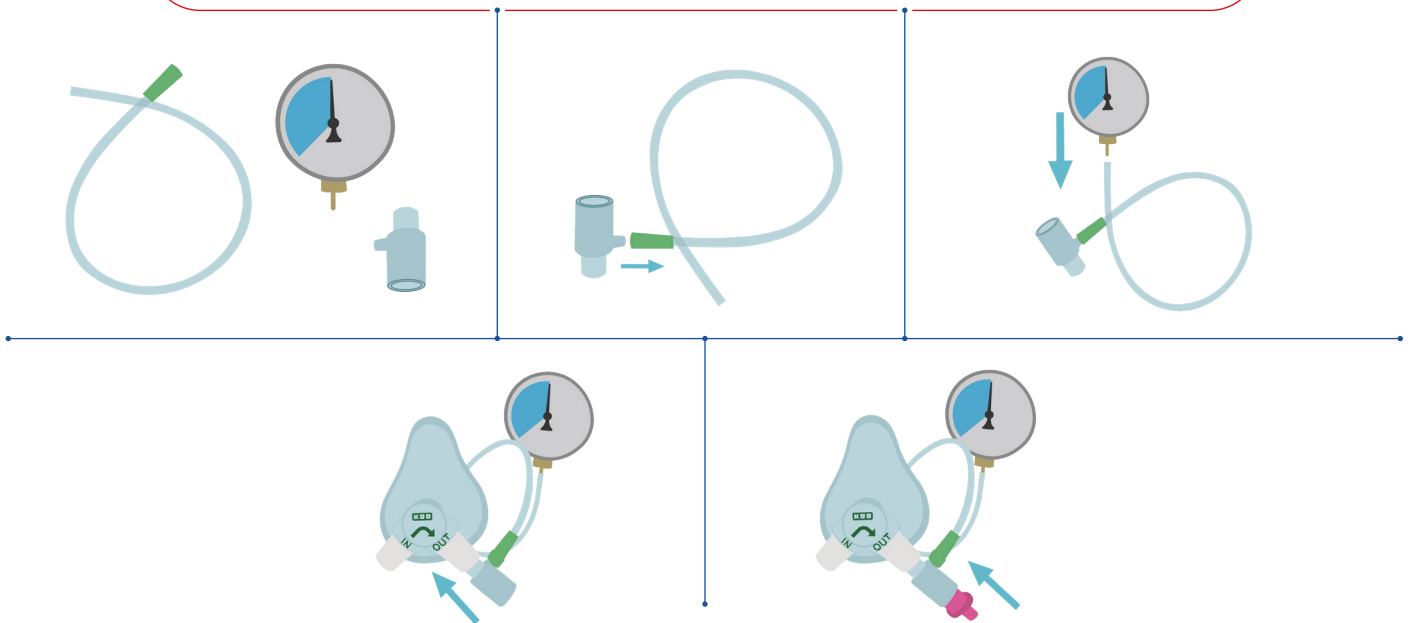

Kit assembly with ECO manometer

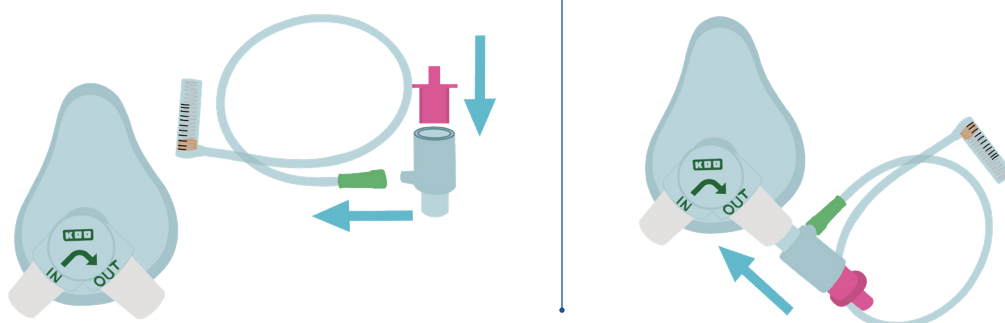

KOO MEDICAL



www.koomedical.com

KOO (SHANGHAI) Industries Co.Ltd. 100 Zhongde Road, Xiaokunshan Town. Songjiang Shanghai 201614 China

Distributed in EMEA Koo by Europe srl Viale delle Industrie 5 - 20044 Arese (MI) - Italy info@kooeurope.com – phone +39 02 93589495

Distributed in Asia Pacific by Koo Asia Ltd. Rm 2113, Level 21, Landmark North 39 Lung Sum Avenue, Sheung Shui, N.T., Hong Kong.
Phone: +852 2110 4780 - Fax: +852 2110 4783